


WP 510.02

Set of 5 torsion specimens, CuZn

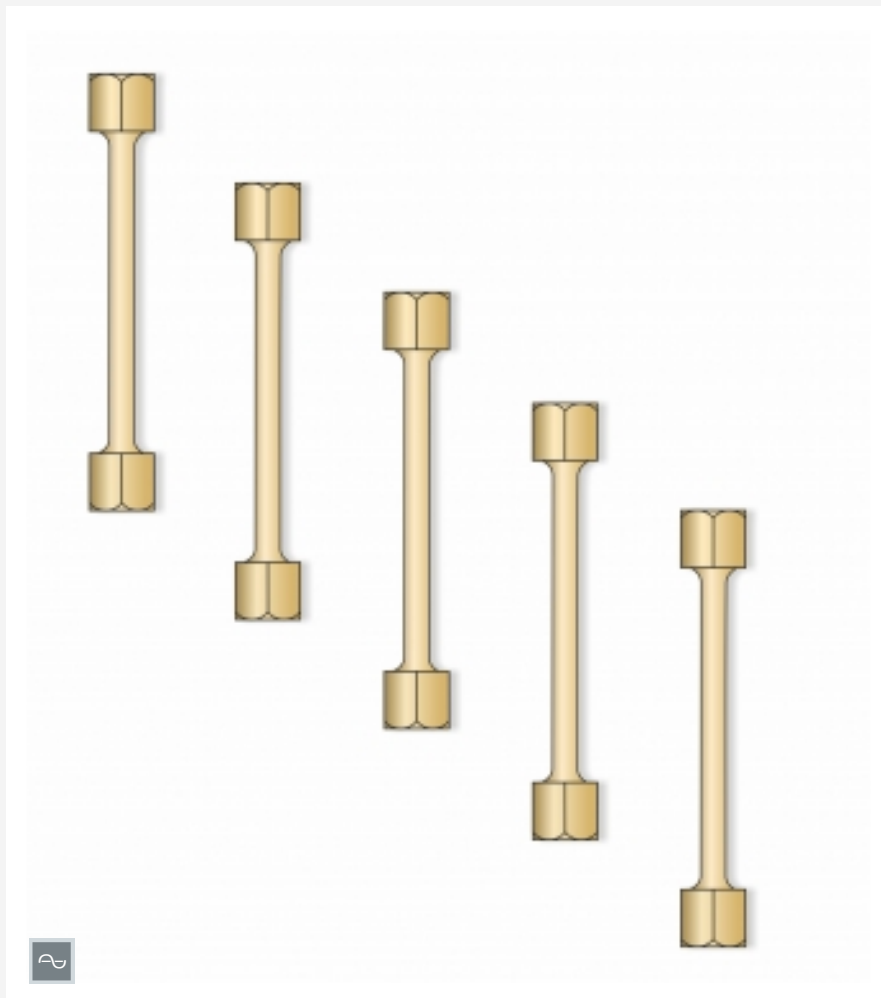

Learning objectives/experiments

- torsion tests on the WP 510 experimental unit

Specification

- [1] torsion specimens for torsion tests
- [2] 5 specimens made of brass
- [3] accessory for WP 510 Torsion test 200Nm, motor drive

Technical data

- Specimens, brass
- length: 146mm
 - clamping length: 100mm
 - specimen diameter in measuring cross-section: 9mm
 - clamp fixing: SW19

Weight: approx. 750g

Scope of delivery

- 1 set of specimens (5 pieces)

Description

■ torsion specimens of round cross-section as accessories for the WP 510 experimental unit

This set of specimens is available as an accessory for the WP 510 experimental unit.

The specimen is installed in the experimental area of WP 510. The ends of the specimens have 19mm hexagons for attachment.

In the experiment, metallic specimens are twisted until they are destroyed by a typical shear fracture. In contrast to the tensile and bending tests, the torsion causes a two-axis stress state. Shear stresses are produced in the specimen which deform it. The deformation is expressed as twisting of the two ends of the specimen towards each other.

The set contains five specimens made of brass.

WP 510.02

Set of 5 torsion specimens, CuZn

Required accessories

020.51000

WP 510

Torsion test 200Nm, motor drive